


**Self-closing doors**  
Doors left open accidentally are a thing of the past. (Except BL 100 PRO-ACTIVE, FROSTER BL 650 PRO-ACTIVE).



**Perspective**  
The optional equipment with a glass door and LED illumination avoids unnecessary door opening.

# BL 520

## PRO-ACTIVE

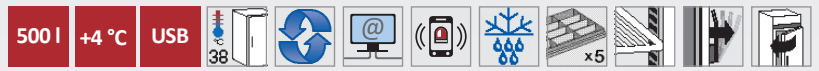

according to DIN 58371

### ■ PRO-ACTIVE- Control:

Permanent, proactive monitoring of the performance data and alerts in the case of deviations – so that you can take countermeasures in good time before a fault threatens your chilled goods; world's most accurate temperature control in refrigerated areas thanks to two standard PT-1000 sensors.

■ **External housing** made from galvanised sheet steel (rust-proof), with white, anti-scratch powder coating. Length of the plug cable: approx. 3 m.

■ **Equipped with castors.**

■ **Interior** made from smooth aluminium with colourless protective coating. Gibs (every 15 mm) for flexible interior arrangement.

■ **Interior consists** of five drawers of aluminium, height adjustable on telescopic runners with stops. Stacking height: 23 cm.

■ **Extra-thick 70-mm insulation**, made from high-quality, compression-moulded and environmentally friendly material. Energy saving.

■ **Self-closing door** with easy-to-replace plastic magnetic seal frame, lockable.

■ **Door stop** on the right-hand side (by default, see illustration), or the left-hand side, can be retrofitted.

■ **Forced-air cooling** with axial blower, switches off automatically when you open the door, ensures a uniform temperature and minimizes temperature deviations.

■ **Automatic fast defrosting** thanks to time-limited and temperature-monitored reversal of the refrigerant circuit.

■ **Condensate evaporates** in stainless steel heated dish beneath the refrigerator floor.

■ **Key switch** protected control panel against tampering. Membrane keyboard with digital temperature display and minimum/maximum temperature memory.

■ **Warning functions** with visual and audible alarm signal in the case of power failure and temperature deviations. In the case of power failure, the monitoring unit remains in operation for approx. 30 hours on battery power. Door open alarm after 60 seconds.

■ **Alarms** can be forwarded using a potential-free contact (e.g., to a mobile phone with optional GSM-Module or to a control centre).

■ **Data documentation** can be read via USB interface with KIRSCH Datanet software.

■ **Antifreeze** against sub-zero temperatures.

■ **Ventilation-enforced refrigerating machine**, vibration free, hermetically sealed, energy saving, low-noise, easy to service. Inlet air and exhaust air flow through the ventilation slits at the front.

■ **Low-noise compressor** reduces noise to 40 dB (corresponds to noise level I; sound perception: whisper).

## Specification

|  |   |
|--|---|
| Capacity   | 500 litres  |
| Temperature setting  | +4 °C   |
| Voltage  | 220 – 240 V, 50/60 Hz   |
| Power consumption  | 250 watts   |
| Normal consumption   | 1.5 kWh/24 h  |
| Admissible ambient temperature                                     | from +10 °C to +38 °C   |
| Heat emission (max.)   | 554 watts   |
| Exterior dimensions (including door handle and distance from wall) | w x d x h = 77 x 76 x 195.5 cm  |
| Interior dimensions  | w x d x h = 62 x 56 x 140 cm (usable width 2 cm, usable depth 11 cm usable height 13 cm less) |
| Exterior dimensions with door open at 90°                          | w x d = 77 x 144 cm   |
| Drawer inner dimensions  | w x d x h = 56 x 39.5 x 10 cm   |
| Max. load drawer   | 50 kg   |
| Weight   | Net weight 138 kg, gross weight 161 kg  |

## Optional equipment

■ **Glass door with lock.**

■ **LED illumination** mounted at the top.

■ **Adjustable feet** to compensate for uneven floors.

■ **Cooling machine** with cold water connection.

■ **Outer housing** in stainless steel 4301, longitudinal brushed.

■ **Additional drawer** made of aluminium.

■ **Drawer insert** (set) with 4 adjustable length and 25 cross dividers per drawer.

■ **Additional length and cross dividers.**

■ **Drawer** made of stainless steel instead of aluminium.

■ **Shelf** 59.7 x 45 cm on rails or overlays.

■ **Baskets** on rails 60 x 45 x 10 cm.

■ **Aluminium tray** with rails 60 x 40 x 2,5 cm.

■ **Condensate container** can be emptied manually if no condensation is required (e.g., in an operating theatre).

■ **Vario Suport**, Storage Shelf on roller runners.

■ **GSM-Module** for alarm text message transmission (e.g., to a mobile phone).

■ **Temperature documentation:**

- PC-KIT-NET

- 4-pin PT 100 or PT 1000 temperature sensor, class 1/3 B, including built-in reference body.

- Disc-type pen recording thermometer, with 7-day measuring cycle, incorporated, with 100 recording discs and 5 felt tip pens.

- Mechanical pen-recording thermometer with two waxed paper strips.